



## *Poz' In' Form®*

---

*“ Proper positioning is crucial in the prevention of pressure ulcers. Due to the material it is made up of, the Poz' In' Form® cushions offer all the required features in bedsore prevention. The combination of a micropore fabric and a bi-elastic and breathable Pharmatex or Lenzing FR® cover offers maximum breathability of the cushion. The filling made of thin microfiber perfectly fits the body shape ensuring optimal pressure distribution and contributes to a significant reduction in friction and shearing forces. ”*

## Effective Positioning

Poz' In' Form® cushions are designed to perform various support and comfort functions for long-term bedridden people. The microfibers prevent from creep effects that are regularly observed in classical positioning and from early deterioration of the support. The various forms, sizes and models of Poz' In' Form® allow for a wide range of positioning possibilities.

## Positioning Comfort

The highly effective positioning cushions reinforce proper positioning of the long-term bedridden client by improving abnormal posture, providing even weight distribution and by preventing pressure sores, while implementing security of use and limiting microbial infections.

## Pharmatex or Lenzing FR® Cover

Poz' In' Form® elements and modules are entirely removable and washable. They are delivered with a Pharmatex or Lenzing FR® cover.

The Pharmatex bi-elastic cover is made of polyurethane and is recommended for use in hospitals. Pharmatex is soft to the touch and very comfortable next to the skin, while still easy to clean.

The Lenzing FR® cover is made of viscose, coated with a waterproof polyurethane film. Lenzing FR® is a major innovation in positioning of the long-term bedridden, preventing maceration. Hygroscopic properties of viscose allow moisture absorption. Moreover elasticity and flexibility of Lenzing FR® are particularly suitable for problems regularly observed in classical positioning.





<b>Code with Pharmatex Cover</b>	7009-0001-269
<b>Code with Lenzing FR® Cover</b>	7009-0001-279



<b>Code with Pharmatex Cover</b>	7009-0551-269
<b>Code with Lenzing FR® Cover</b>	7009-0551-279
<b>Length</b>	56 cm
<b>Depth</b>	32 cm
<b>Height</b>	28 cm



<b>Code with Pharmatex Cover</b>	7009-1251-269	7009-1252-269	7009-1253-269
<b>Code with Lenzing FR® Cover</b>	7009-1251-279	7009-1252-279	7009-1253-279
<b>Length</b>	35 cm	55 cm	40 cm
<b>Depth</b>	25 cm	40 cm	15 cm

## Poz' In' Form Kit

- Transport bag
- Universal Cushion 55cm x 40cm
- Half moon cushion
- Semi-lateral decubitus cushion
- Cylindrical cushion
- Ring cushion
- Heel cushion with strap
- Occiput head cushion
- Pharmatex covers or Lenzing covers



### Effective Positioning

Poz' In' Form® cushions are designed to perform various support and comfort functions for long-term bedridden people. The microfibers prevent from creep effects that are regularly observed in classical positioning and from early deterioration of the support. The various forms, sizes and models of Poz' In' Form® allow for a wide range of positioning possibilities.





## Poz' In' Form Delta

- This cushion adapts to all morphologies. Placed under knees, it supports the knees' joint, prevents the patient from sliding downwards and reduces friction and shearing effects on the sacrum's area.
- It can be used in a semi-fowler position and helps to maintain a perfect position while relieving pressure in sensitive areas (sacrum, ischial tuberosities, heels) and shifting it to low-risk zones (calves, thighs and back).

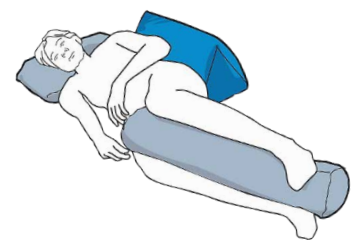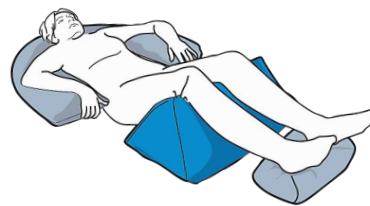


## Poz' In' Form Universal

- Thanks to the three sizes, the universal cushion can be used in many different ways. In lying down the 55 x 40 cm size can be used for the knee flexion and a moderate hip abduction.
- The smaller cushions help to hold the patient in optimal positioning thus reducing bedsores risks and abnormal postures.

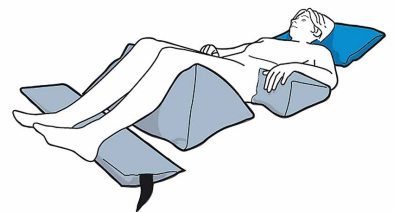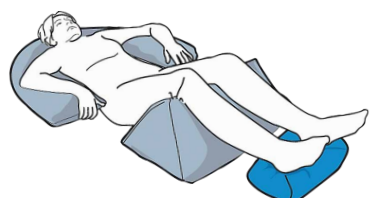

## Poz' In' Form Cylindrical


- Placed between legs, it reduces friction and shearing effects on thighs and knees' area.
- It can also be used under ankles to relieve pressure.

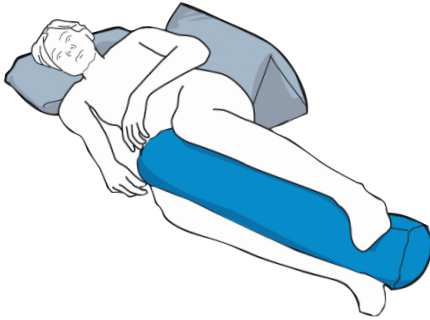

<b>Code with Pharmatex Cover</b>	7009-1501-269
<b>Code with Lenzing FR® Cover</b>	7009-1501-279
<b>Length</b>	75 cm
<b>Diameter</b>	18 cm

## Poz' In' Form Abduction


- It is used in crotch area when the patient is lying down to hold hips abducted and legs apart to reduce friction and shearing effects between legs, on the inner sides of knees and calves.

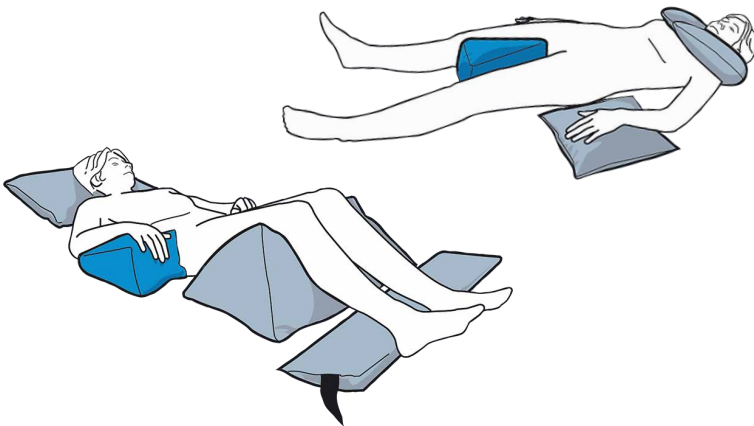

<b>Code with Pharmatex Cover</b>	7009-1551-269
<b>Code with Lenzing FR® Cover</b>	7009-1551-279
<b>Length</b>	35 cm
<b>Small Width</b>	10 cm
<b>Large Width</b>	30 cm
<b>Height</b>	15 cm

## Poz' In' Form Half-Moon

- In supine or semi-fowler position it supports and holds the back in proper position. It can also be used in (wheel) chairs for a comfortable and upright sitting position.

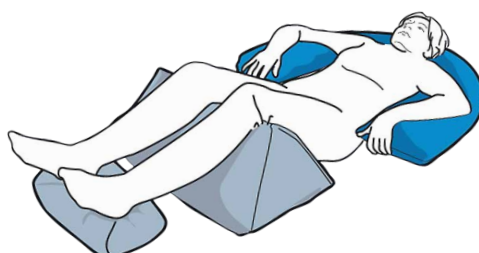

<b>Code with Pharmatex Cover</b>	7009-1851-269
<b>Code with Lenzing FR® Cover</b>	7009-1851-279
<b>Length</b>	200 cm
<b>Width</b>	30 cm



<b>Code with Pharmatex Cover</b>	7009-3001-269
<b>Code with Lenzing FR® Cover</b>	7009-3001-279
<b>Length</b>	185 cm
<b>Width</b>	35 cm




	Half Ring	Ring
<b>Code with Pharmatex Cover</b>	7009-4051-269	7009-4081-269
<b>Code with Lenzing FR® Cover</b>	7009-4051-279	7009-4081-279
<b>Length</b>	105 cm	135 cm
<b>Width</b>	22 cm	20 cm



<b>Code with Pharmatex Cover</b>	7009-4751-269
<b>Code with Lenzing FR® Cover</b>	7009-4751-279
<b>Length</b>	50 cm
<b>Width</b>	40 cm
<b>Height</b>	9 cm


## Poz' In' Form Semi-Lateral Decubitus

- A multi-purpose cushion, it comfortably supports shoulders, back, hips and knees when the patient is lying in lateral position. Moreover, it reduces friction and shearing effects.

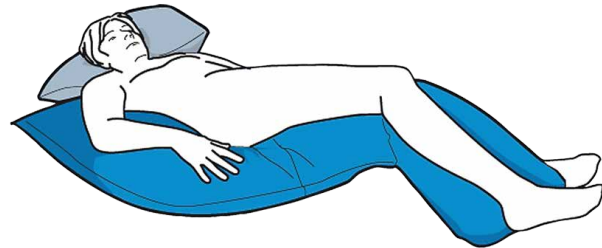


## Poz' In' Form Ring & Half Ring

- In combination with the cylindrical cushion, the half ring cushion reduces friction and shearing effects between legs, on the inner sides of knees and thighs and holds hips abducted. It can also be used under the neck.
- This cushion can be used in addition to universal, triangle or cylindrical cushions to support the posture correction. Put around the neck, this cushion comfortably supports the patient's head.

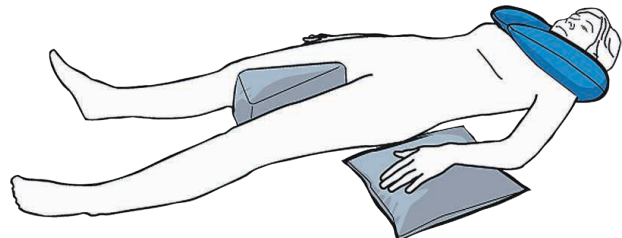


## Poz' In' Form Occiput

- The hollow ensures an optimal positioning and a pressure relief on the occipital area. The combination of memory foam and microfiber allows a regular relief on the contact zones.

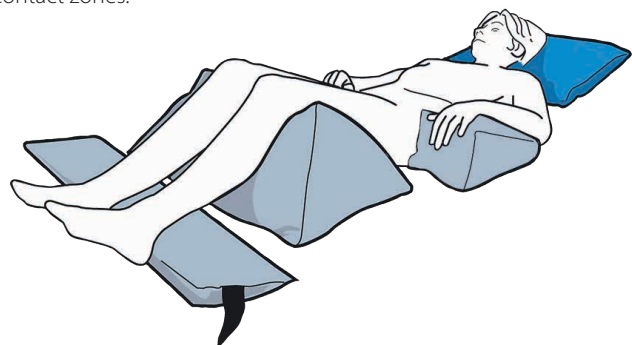

## Poz' In' Form Heel Cushion


- It is placed at the bottom of the bed and fixed with an adjustable strap.
- The pressure on the heels is completely discharged.
- The heel cushion holds the ankles in proper position.

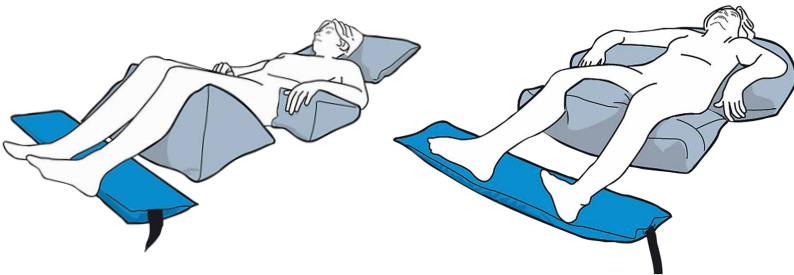

<b>Code with Pharmatex Cover</b>	7009-5201-269
<b>Code with Lenzing FR® Cover</b>	7009-5201-279
<b>Length</b>	85 cm
<b>Diameter</b>	25 cm

## Poz' In' Form Heel Wedge


- This cushion contributes to shift the pressure points from the high-risk heel area to low-risk zones (calves and knees).
- The heel is completely off the mattress and the weight of the lower leg is distributed evenly across the cushion. It holds the ankles in proper position and leads to a better blood circulation.
- Its size allows pressure relief on the Achilles' heel.

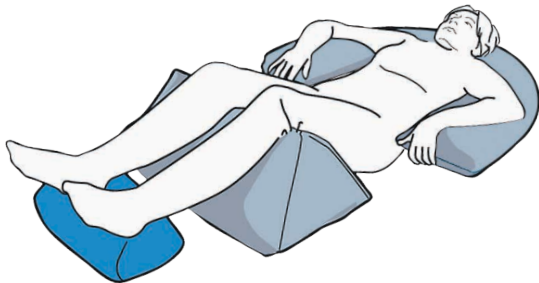

<b>Code with Pharmatex Cover</b>	7009-5401-269
<b>Code with Lenzing FR® Cover</b>	7009-5401-279
<b>Length</b>	48 cm
<b>Width</b>	30 cm
<b>Height</b>	12 cm

## Poz' In' Form Hand Cushion

- This positioning cushion for the hand reduces the pressure on the skin on the zones with risk of pressure sore. It is advisable after surgical operations.
- It reduces the number of pressure sensitive bearing points and allows a discharge of the elbows area. Its size allows pressure relief on the wrist.

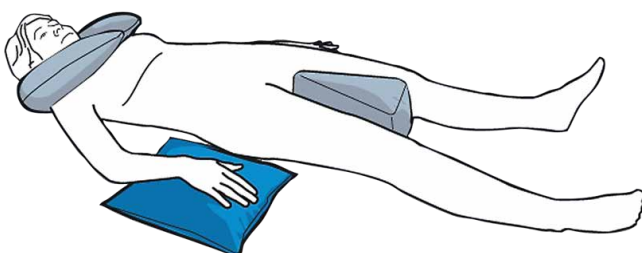

<b>Code with Pharmatex Cover</b>	7009-5411-269
<b>Code with Lenzing FR® Cover</b>	7009-5411-279
<b>Length</b>	45 cm
<b>Width</b>	23 cm
<b>Height</b>	8.5 cm





*Shaping better lives.®*

+64 3 307 9790  
0800 543 343  
[solutions@medifab.com](mailto:solutions@medifab.com)